

# EGIS

**EGIS**


Esophageal  
Stent

## EGIS **Esophageal** Stent

E-PTFE Membrane Stent to Esophageal



Full Anti-Reflux Silicone

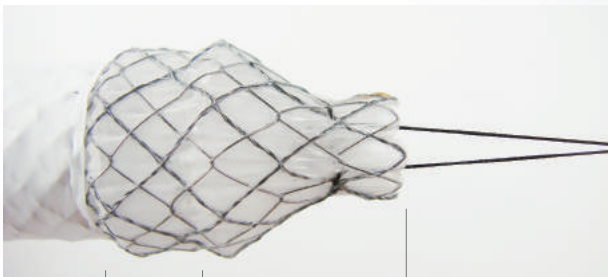

Body

Head  
-1

Head  
-2

# EGIS

## Esophageal Stent

# EGIS Esophageal Stent

E-PTFE Membrane Stent to Esophageal

### Esophageal Stent



EGIS Esophageal Stent is unique in design with anti-reflux option. The teflon coating is used as it lasts longer than silicone when in contact with stomach acid. Double-stepped shoulders ("DUMBBELL" shaped) minimize stent migration. The e-PTFE membrane, usually used as an artificial vessel in stent graft, provides a firm but flexible feature which will not perforate when it is implanted or removed.



### Specification

Delivery System 15/18Fr, 70cm

Ordering No. First two\*\* : Diameter(mm) / Last two\*\* : Length(cm)

TYPE	Ordering No.	Diameter(mm)		Length(cm)	
		Body/Head-1/Head-2	body	Total	
Bare	EBA****07				
Full Covered Silicone	ESC****07				
Full Covered e-PTFE	EEC****07	16/20/24	4	9	
Full Anti-Reflux Silicone	ESA****07	18/22/26	6	11	
Full Anti-Reflux e-PTFE	EEA****07	20/24/28	8	13	
Double Covered Silicone	EDC****18		10	15	
Partial Covered Silicone	ESP****07		12	17	
Partial Covered e-PTFE	EEP****07	18/22/26			
		20/24/28			