

Ascendia MedTech AB introduces  
a brand **new idea** in cover technology:

# Round Seams

No trapped air

No weld ears

This is  
TopFit®

No seams that  
hinder imaging

This is  
an ordinary  
cover

## TopFit® Sterile Ultrasound Covers

### The TopFit® Advantages

- An ultrasound cover with a circular shaped seam - patent pending.
- Less image distortion.
- No ears - less risk for penetrating the cover with the biopsy needle.
- Tailor made cover with snug fit, and a minimal seam minimizes air bubbles and risk of leakage.
- Made of soft and pliable latex free material.



# TopFit® Sterile Ultrasound Covers

## Think outside the box



The box end ultrasound cover, also known as the three-dimensional tip 3D, has become the standard for non-latex ultrasound covers. Even though most commonly used transducers are not really box shaped...

One drawback with the box end is its mechanical properties, leading to a loose fit over the transducer, allowing for air bubbles to degrade the image quality. The traditional placing of seams also creates weld ears that can cause image distortion and increase the risk of penetrating the cover with the biopsy needle.

The Ascendia research team has now created the logical and innovative solution: TopFit®, the cover with a circular shaped seam. Improving imaging properties and ease of use; there is no need for alignments since the transducer, due to the roundness of the cover, aligns automatically.

### Ascendia MedTech AB. Better imaging properties by better design.

Ascendia MedTech AB has since the start in 1985 been serving the ultrasound community with products that use innovative design and materials in order to improve ease of use and image quality. One good example is that our covers ensure a better fit, a safer and more accurate procedure and good value. Because that is what you need.

### Cover Shapes & Transducer Types

#### Mid/Large Size Linear/Curved TXD


Article number: 105-7324 Width: 10,2 cm Length 61 cm  
Article number: 106-7324 Width: 10,2 cm Length 91,5 cm  
Article number: 107-7324 Width: 10,2 cm Length 122 cm  
Article number: 108-7324 Width: 10,2 cm Length 147 cm  
Article number: 109-7324 Width: 10,2 cm Length 244 cm  
Article number: 105-7424 Width: 8,9 cm Length 61 cm  
Article number: 106-7424 Width: 8,9 cm Length 91,5 cm



The circular shaped seam also makes this cover a perfect fit for most commonly used ultrasound probes on the market.

#### Mid Size Curved TXD

Article number: 105-7524 Width: 14 cm Length 61 cm  
Article number: 106-7524 Width: 14 cm Length 91,5 cm  
Article number: 107-7524 Width: 14 cm Length 122 cm  
Article number: 108-7524 Width: 14 cm Length 147 cm  
Article number: 109-7524 Width: 14 cm Length 244 cm

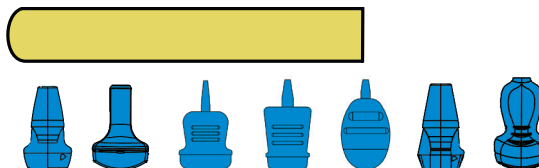

### Kit contents

24 Procedure kits per box, each containing:  
1Ea TopFit® Cover  
1Ea SafeScan™ Sterile Gel  
2Ea fasteners, non-latex  
2Ea rubber bands, non-latex

Article number 109-7324 and 109-7524 contains  
12 Procedure kits per box.

Distributed by:



ASCENDIA MEDTECH AB  
Svärdvägen 21 • SE-18233 Danderyd • Sweden  
Phone + 46-8-444 5400  
Email: info@amedic.se  
www.amedic.se