


SUBINTIMAL ANCHOR TECHNIQUE



SIDE BRANCH PROTECTION



SIDE BRANCH ACCESS

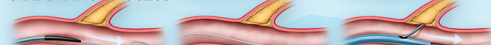


SIDE BRANCH ACCESS THROUGH STENT STRUT




SIDE BRANCH ANCHOR TECHNIQUE



RECROSS 360° ACCESS

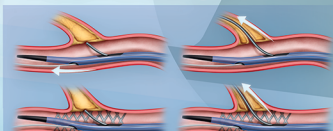


NHANCER Rx

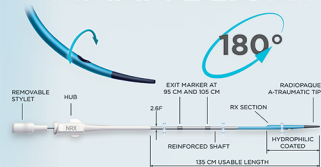

DUAL LUMEN PROCEDURAL TECHNIQUES



RECROSS

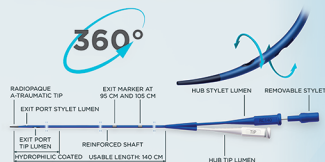


REVERSE WIRE TECHNIQUE




SELECTIVE CONTRAST INJECTION




ANTEGRADE DUAL ACCESS (ADA)

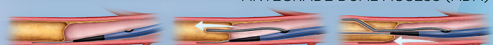

SUBINTIMAL GUIDEWIRE REDIRECTION (SGR)

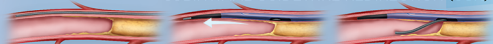

HEMATOMA DECOMPRESSION

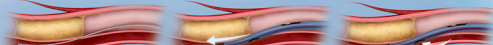

IMDS

INTERVENTIONAL MEDICAL
DEVICE SOLUTIONS



NHANCER Rx



ADVANCING
DUAL GUIDEWIRE
PROCEDURES

WWW.IMDS.NL

REVERSE WIRE

SIDE BRANCH ACCESS
THROUGH STENT STRUT

SUBINTIMAL ANCHOR

SELECTIVE CONTRAST INJECTION

SIDE BRANCH
ANCHOR

NHANCER Rx

DISTAL SECTION

FULL-LENGTH VARIABLE BRAID


SOFT TAPERED TIP




DIMENSIONAL DETAILS




POSITIONING AND TORQUE CONTROL

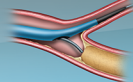

TRUE VISUALISATION
ON BOTH EXIT PORTS

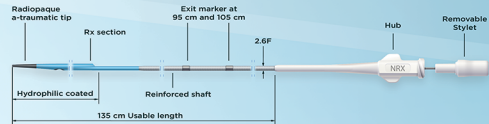

OVAL DISTAL SHAFT MINIMIZING
CROSSING PROFILE & SUPERIOR PUSH
EFFICIENCY FOR SUPPORT
OF CORONARY ACCESS

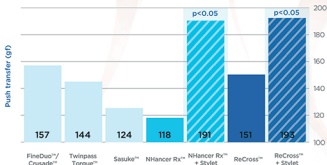

RECROSS

6F GUIDING CATHETER COMPATIBLE




Distal crossing surface (mm²)



NHANCER Rx

LESION CROSSING
SUPPORTED BY BEST IN CLASS
TIP PENETRATION

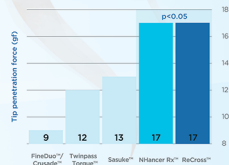

SCAN CODE FOR
INSTRUCTION VIDEOS



GUIDEWIRE COMPATIBILITY OTW LUMEN	0.014"	DISTAL SHAFT COATING	HYDROPHILIC (NDURANCE)
GUIDEWIRE COMPATIBILITY RX LUMEN	0.014"	RX LUMEN LENGTH	18 cm
SHAFT INNER DIAMETER	0.48 mm	USABLE LENGTH	135 cm

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.
Caution: Federal Law (U.S.) restricts this device to sale by or on the order of a physician.
In-vitro test data on file at IMDS.
Bench test data may not be indicative of clinical performance.
NHancer™ is a trademark of IMDS.

© IMDS 2023

ITM200207 R6.0



Ceintuurbaan Noord 150
9301 NZ Roden
The Netherlands
sales@imds.nl
+31(0)560 8200230
Phone
Fax +31(0)560 8200231